



Patient Guide

First Medical Source (FMS) Pain Management System

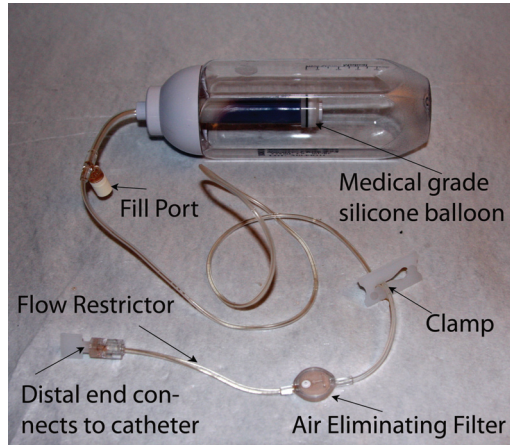


CAUTION: FMS Pain Management system consists of the pump, the catheter(s), and all the necessary instructions. Use the system only as directed by your healthcare provider.

TO START MEDICATION DELIVERY: Open the clamp

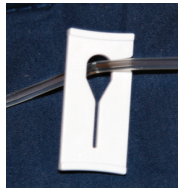

OPEN CLAMP

TO STOP MEDICATION DELIVERY: Close the clamp

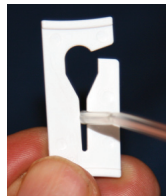

CLOSED CLAMP

HOW THE PUMP PERFORMS: Under elastic pressure, the balloon inside the hard case pump will deliver medicine to the surgical site. For the delivery rate to be accurate, the flow restrictor at the end of the tubing (see above picture) must be in contact with your skin. Tape it to your skin if needed.

DURING MEDICATION DELIVERY WATCH FOR THE FOLLOWING:

- The fluid level in the pump is reduced
- The system is not leaking
- No air has been trapped in the medication reservoir
- The tubing is not kinked
- System connections are all secure

If you observe any problem, immediately contact your healthcare provider.



You may continue your active life during the infusion by wearing the pump. You may either clip the pump to your belt, or place it in the carry pouch (if provided).



As the pump delivers your medication, the balloon inside the reservoir will slowly reduce size. You may monitor the medication delivery by observing the fluid level using the level indicator lines on the pump.

MEDICATION DELIVERY IS COMPLETE

If your healthcare provider has instructed you to remove the catheter(s) after the delivery is completed, follow their instructions, keeping in mind the following points:

TO REMOVE THE CATHETER

1. Gently remove the dressing



2. Hold the catheter near the skin and carefully pull it out. It should not be painful or hard to remove. If it is stuck, or difficult to remove, STOP and call your healthcare provider. DO NOT cut the catheter. DO NOT try to remove the catheter forcefully.



3. After you removed the catheter, place a bandage over the catheter site. Make sure you can see the black mark at the end of the catheter.



4. If you do not see a black mark at the end of the catheter (this may mean you have left a piece of the catheter in place) or if you see any redness or feel warmth near the incision site, call immediately your healthcare provider.

5. If you have more than one catheter, repeat these steps for each catheter.



IMPORTANT INFORMATION

MEDICATION

Name:

Dosage:

Infusion time:

**PROVIDER
MEDICATION**

Name:

Doctor's name

Phone number:

Special Instructions:

Questions to ask or concerns to
share with your doctor:



First Medical Source

Reducing Cost of Healthcare Worldwide

Patient Guide

First Medical Source (FMS) Pain Management System



PRECAUTIONS:

- Do not reuse
- Do not immerse in water (shower, bath or swim)

Report any of the following symptoms to your provider:

- Increased pain
- Redness, swelling, or discharge in the incision area
- Numbness, dizziness, or drowsiness
- Nausea
- Metallic taste in your mouth
- Other Side-Effects

To learn more visit our website:

www.firstmedicalsourc.com

CAUTION: Federal law restricts this device to sale by or on the order of a physician

Manufactured for: First Medical Source, LLC, Laguna Niguel, CA 92677

© 2010 By First Medical Source, LLC, All rights reserved.